Piezon Master 400

Operating Instructions

Mode d'emploi

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English

Dear customers,

These operating instructions include four languages: English, French, Italian and Spanish. The English instructions are on pages E-1 to E-19.

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Technical modifications

EMS reserves the right to modify the technique, accessories, operating instructions or contents of the original Piezon Master 400 set due to technical or scientific improvements.

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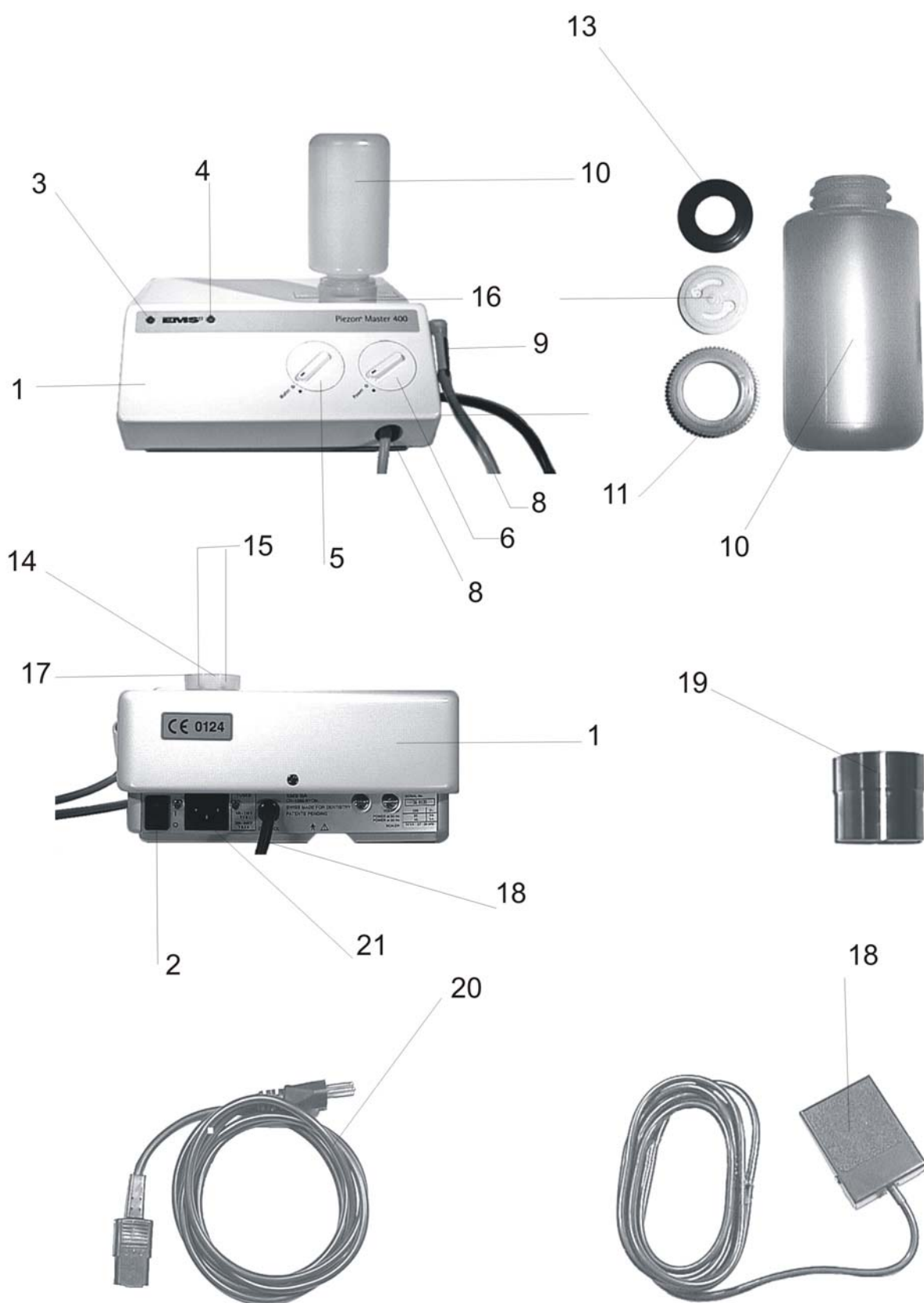
Piezon Master 400

Dear customers,

Thank you for buying a new EMS unit. Your Piezon Master 400 meets the highest quality and safety standards..

This ultrasonic unit is for use in dentistry with EMS ultrasonic instruments.

EMS supplies ultrasonic instruments for the following indications:

- Removal of supragingival calculus
- Removal of subgingival calculus
- Periodontal treatment
- Preparation of approximal cavities
- Luting tooth-shaded inlays and onlays with highly thixotropic, dual-curing cements
- Removal of crowns, bridges, inlays and posts as well as condensing gutta-percha
- Preparation, cleaning and irrigation of root canals
- Retrograde preparation of root canals

For further details on the various ranges of indications and handling the instruments, please refer to the handling instructions which are supplied with your Piezon[®] systems or are available from the distributor in your country.

Please read this before beginning work!

These operating instructions are intended to help you install, connect and operate your Piezon Master 400 properly.

Please read these operating instructions carefully - they explain all the most important details and procedures. Please pay special attention to the safety precautions.

Always keep these instructions close at hand.

To prevent injury to people and damage to property, please heed the relevant warnings and remarks. They are marked as follows:

- Danger
Where people are at risk
- Caution
Where property may be damaged and the environment harmed
- Please note
Useful additional information and hints

We would be pleased to answer your questions or listen to your suggestions about the Piezon Master 400. For technical support, please call our EMS service staff or your dealer. The telephone numbers and addresses of the EMS facility in your area are on the last page of these operating instructions.

We wish you lots of success!

EMS

Contents

Safety precautions	E-4
Controls and components of the Piezon Master 400	E-5
Connecting the Piezon Master 400	E-5
Operating the Piezon Master 400	E-6
"Water" and "Power" knobs	E-6
Rocker switch - ENDO or IRRIGATION/SCALING - and footswitch	E-6
Handpiece holder	E-6
Filling and fitting the bottle	E-7
Torque wrench	E-8
Disinfecting, cleaning and Sterilizing the Piezon Master 400	E-8
Regular care and maintenance	E-8
Piezon Master 400 accessories	E-9
EMS Service	E-9
Storing the Piezon Master 400 when out of use for extended periods	E-10
Disposing of the Piezon Master 400, accessories and packaging	E-10
Technical data	E-11
Symbols	E-12
Electromagnetic compatibility	E-13
Electromagnetic emissions	E-13
Electromagnetic immunity	E-14
Electromagnetic immunity, RF portable equipment	E-15
Recommended separation distances	E-16
Cables length	E-17
Troubleshooting	E-18

Safety precautions

EMS accepts no liability for direct or consequential injury or damage resulting from improper use, arising in particular through the non-observance of the operating instructions, or improper preparation and maintenance. The provision of technical documents does not imply any authorization for any adjustment, repair or modification of this EMS product.

Danger - New or repaired handpieces or instruments must be sterilized before first use
Before first use, clean, disinfect and sterilize new or repaired handpieces and instruments which you wish to use with the Piezon Master 400. Please refer to the information provided in the operating instructions.

Danger - Use for the intended purpose only
Failure to observe the operating instructions may result in the patient or user suffering serious injury or the product being damaged, possibly beyond repair. Before using the product, make sure that you have studied and understood the operating instructions. This also applies to any product used with the Piezon Master 400.

Danger - For use by qualified and trained personnel only
The Piezon Master 400 must not be used by other than properly qualified and trained personnel, in medical facilities.

Caution - Always examine the unit and accessories for damage before commencing treatment
Always examine the unit and accessories for damage before commencing treatment. A damaged Piezon Master 400 or damaged accessories must not be used and must be replaced.

Danger - Do not install where there is a risk of explosion
The Piezon Master 400 is not intended for operation in presence of inflammable anaesthetics or gases.

Danger - Contraindications
Ultrasonic oscillations may prevent cardiac pacemakers functioning properly. Therefore, EMS recommend that patients with a cardiac pacemaker should not be treated with the Piezon Master 400.

Caution - Check the safety status after repair
The safety status of the Piezon Master 400 and accessories must comply with the applicable standard for this type of product after it has undergone any type of repair or modification. A report must be drawn up in respect of every inspection, repair or modification.

Caution - This product needs special precautions regarding electromagnetic compatibility and needs to be installed and put into service according to the section "Electromagnetic compatibility" provided in this document.

Caution - Certain types of mobile telecommunication equipment could potentially interfere with this product. The separation distances recommended in the section "Electromagnetic compatibility" of this document must be taken into account

Danger - This product should not be used adjacent to or stacked with other equipment. If adjacent or stacked use is necessary, this product should be observed to verify normal operation in the configuration in which it will be used.

Danger - The use of accessories and cables other than those specified or sold by EMS as replacement parts, may result in increased emissions or decreased immunity of this product.

Controls and components of the Piezon Master 400

- 1 Housing
- 2 Power switch
- 3 Equipment-on lamp, green
- 4 Control lamp, yellow
- 5 "Water" knob for controlling the water flow rate
- 6 "Power" knob for controlling the ultrasonic power
- 7 ENDO or IRRIGATION/SCALING rocker switch
- 8 Handpiece hose with quick-link coupling for Piezon handpieces
- 9 Handpiece holder
- 10 Irrigating liquid bottle
- 11 Screw-cap for bottle
- 12 Bayonet face plate for the screw-cap
- 13 Red washer for the screw-cap
- 14 Bottle holder
- 15 O-ring for the small nipple in the bottle holder
- 16 Silicone washer for the bottle holder
- 17 Bottle holder cap
- 18 Footswitch
- 19 Torque wrench
- 20 Power cable
- 21 Power socket

Connecting the Piezon Master 400

Please note - It is advisable to keep the original packaging of the unit. Keep the packaging for storing or transporting the unit.

You require a power supply in order to connect the Piezon Master 400.

Use the power cable (20) to connect the unit to the power supply. Before connecting, check that the rated voltage specified on the back of the unit corresponds with the local line voltage of your dental office. The rated voltage of the unit must be the same as the local line voltage, otherwise the unit can be damaged. Plug the power cable (20) into the power socket (21) on the back of the unit and then into a mains power socket.

Operating the Piezon Master 400

Switch the Piezon Master 400 "ON" with power switch (2). The green equipment-on lamp (3) on the front of the unit will illuminate.

Please note - Read the operating instructions for the appropriate EMS system. They cover the "Water" knob (5) and "Power" knob (6) settings for every instrument in detail.

"Water" and "Power" knobs

The "Water" knob (5) and "Power" knob (6) are on the front of the Piezon Master 400. The irrigating liquid flow rate and ultrasonic setting can be adjusted progressively with the knobs (5, 6) from the lowest setting, via the low and medium settings, up to the highest setting. Set the knobs (5, 6) as described in the operating instructions for the system in use.

ENDO or IRRIGATION/SCALING rocker switch and footswitch

The ENDO or IRRIGATION/SCALING rocker switch (7) is on the right hand side of the Piezon Master 400, beneath the handpiece holder (9). The rocker switch has two settings which control the two-stage footswitch (18).

When carrying out treatment with the Piezon Master 400, set the rocker switch to IRRIGATION/SCALING. When the rocker switch is set to this position, the footswitch should be pressed lightly to "irrigate only". Press the footswitch firmly to actuate the "ultrasonics and irrigation".

The "ENDO" rocker switch setting is required for orthograde endodontic treatment only. When the rocker switch is set to "ENDO", the footswitch (18) should be pressed lightly to actuate "ultrasonics and irrigation". Press the footswitch firmly to actuate "ultrasonics only". When the footswitch is in this setting, the yellow control lamp (4) illuminates and the sound of the pump subsides. Please note that in this setting, the irrigating liquid is heated immediately - only use the second footswitch setting briefly. The brief duration of the treatment prevents possible injury to the periodontium.

Handpiece holder

The handpiece holder (9) includes a safety mechanism which prevents the handpiece being operated inadvertently as long as it is in the handpiece holder (9).

Filling and fitting the bottle

The following irrigating liquids can be used, depending on the type of system and treatment:

- o physiological saline solution
- o distilled water
- o sodium hypochlorite 0,5%

Caution - The irrigating liquid bottle is not suitable for storing medical solutions.

Danger - Due to sterilization and certain chemicals such as sodium hypochlorite, the irrigating liquid bottle will become brittle and cracked in time. You must check the bottle and screw-cap regularly for signs of damage. A damaged bottle may burst due to the pressure built up in it during treatment. The bottle must be replaced at the first signs of cracking or brittleness.

The irrigating liquid bottle comprises two parts - the bottle and the screw-cap (11). The screw-cap consists of three parts - the threaded ring, the bayonet face plate (12) and the red washer (13). The bayonet face plate (12) retains the bottle in the holder. The red washer (13) should be placed flat in the screw-cap (11). It seals the screw-cap between the bottle and bayonet face plate.

To fill the bottle (10) with irrigating liquid, unscrew the cap. Before closing the bottle, ensure that the red washer (13) has been placed flat in the screw-cap (11). Tighten the cap, but without using excessive force.

Caution - When positioning and removing the bottle, the "Water" knob (5) must be at its lowest setting and the handpiece hose in the handpiece holder.

Remove the bottle holder cap (17) on the unit. Check that the transparent silicone washer (16) is in the bottle holder. Place the bottle (10), with its screw-cap (11) downwards, in the bottle holder (14). The bayonet catch ensures that the bottle will not drip and is airtight. Turn the bottle anti-clockwise until it locks into the bayonet face plate (12) properly. Then turn the bottle (10) a quarter-turn clockwise without using force. This opens the bayonet catches of the bayonet face plate. The bottle is now ready for use.

To replace the bottle, please release the footswitch (18) and place the handpiece hose (8) in the handpiece holder (9). Turn the bottle (10) a quarter turn anti-clockwise and remove it from the bottle holder (14). Pressure builds up in the bottle (10) during treatment. The pressure must be released as, otherwise, the bottle may be damaged. To release the pressure, turn the screw-cap of the bottle (11) a quarter turn anti-clockwise. The air will escape. Tighten the screw-cap again by turning it clockwise, but without using excessive force.

Torque wrench

The torque wrench is used to fix the EMS instruments onto the appropriate handpiece using the correct amount of force. When fixing file holders into place, please use the special wrenches included with those systems. The torque wrench must not be used for this purpose. Please also refer to the operating instructions of the particular systems.

First tighten the instrument gently by hand. Then use the torque wrench (19). Continue turning until you feel resistance on the torque wrench (19). Then turn the wrench another quarter turn. The instrument will then have been tightened optimally.

Disinfecting, cleaning and sterilizing the Piezon Master 400

After treatment, always rinse the Piezon Master 400, with coupled handpiece, in lukewarm water. This prevents the various disinfectants reacting in the sections which they pass through, which would clog and prevent proper functioning. Disinfectants may attack the material which the unit is made of or the handpiece hose if they are exposed to them for longer periods.

To flush through the sections which the liquid passes through, half fill the bottle (10) with water. Adjust the "Water" knob (5) to its highest setting and rinse for more than 20 seconds.

Then pull the handpiece off the handpiece hose (8). Shake the water out of the handpiece and hose. This prevents short circuiting in the handpiece hose plug.

Clean and disinfect the torque wrench (19) after unscrewing it from the handpiece. The torque wrench can be autoclaved at up to 135°C for at least 3 mins.

Regular care and maintenance

The irrigating liquid bottle can be sterilized. Before sterilizing the bottle, clean it with running water. Do not use scouring powder, sponges which scratch, or brushes.

Caution - Do not sterilize the bottle in a chemical vapour sterilizer or hot-air sterilizer as this would damage the bottle.

The bottle (10) can be autoclaved at a maximum of 121°C and 1,0 bar for a maximum of 15 minutes. Unscrew the cap (11) from the bottle (10). Otherwise, excessive pressure will build up in the bottle, which damages it during sterilization. To sterilize the bottle, put it on a tray insert. Do not lay the bottle (10) or screw-cap (11) on bare metal. Heat may be transferred and melt the bottle (10) and screw-cap (11). Check the condition of the bottle and screw-cap after sterilization.

The housing (1) of the Piezon Master should be disinfected with commercially available disinfectant. Do not use scouring powder as it scratches the surface of the housing.

Danger - Please note that the housing of the Piezon Master 400 is only protected against splashed water and must not be sterilized.

If the irrigant liquid bottle (10) can no longer be closed tightly, replace the red washer (13) on the screw-cap.

The bottle holder (14) has two small nipples which engage the openings in the bayonet face plate (12). A small black O-ring on each of these nipples forms a seal between the bayonet face plate and nipple. If the O-rings are damaged or worn, irrigating liquid may leak out. Should this occur, check both O-rings (15) and replace them if necessary. Should irrigating liquid continue to leak out, the bayonet face plate (12) is defective and must be replaced.

Piezon Master 400 accessories

EMS systems are available from dental dealers. Please contact your dental dealer.

EMS-Service

Should your Piezon Master 400 need additional servicing or repairs, please send it to your dealer or our EMS service staff. The telephone numbers and addresses of the EMS service staff covering your area are on the last page of these operating instructions.

The EMS service staff are committed to servicing and repairing your equipment as quickly as possible.

Only have your Piezon Master 400 repaired by EMS or someone authorized by EMS. Only use original EMS spare parts and accessories. Should the unit be opened without authorization or damaged due to non-adherence to the operating instructions, EMS accepts no liability whatsoever. This will also void the guarantee.

Please note - It is best to ship your Piezon Master 400 in the original packaging. It protects your Piezon Master 400 against damage during shipping. Before dispatching your Piezon Master 400, please clean and disinfect it and sterilize the handpiece and instruments as described in the instructions. Always ship your Piezon Master 400 with the complete accessories. When sending your Piezon Master 400 directly to the EMS service department, please include the name and address of your dental dealer. This simplifies processing for us.

Storing the Piezon Master 400 when out of use for extended periods

Should you wish to take your Piezon Master 400 out of use for an extended period of time, please observe the following recommendations:

- fill an irrigant liquid bottle with lukewarm water and rinse the sections where the liquid passes through,
- continue rinsing until no more water runs out of the Piezon Master 400,
- disconnect the unit from the power supply,
- disinfect and clean the housing and handpiece hose,
- pack the unit and all accessories in the original packaging.

This will ensure that, even after being stored for some time, your Piezon Master 400 will still function properly.

Please refer to the storage and transport conditions in the "Technical data" section on page E-11.

Disposing of the Piezon Master 400, accessories and packaging

Keep the original packaging until the Piezon Master 400 is to be disposed of permanently. You can use it for shipping or storing your Piezon Master 400 at any time.

Should you wish to dispose of your Piezon Master 400, please comply with the regulations for disposing of the Piezon Master 400, accessories and packaging which apply in your country.

Technical data

Manufacturer:	EMS SA, CH-1260 Nyon, Switzerland
Model:	Piezon Master 400
Classification:	EN 60601-1: - Class I - Type B - IP 20 93/42 EU: - Class IIa CSA Std. N°125 (*applies to 100/110 V AC units only) - risk class 2G
Mode:	Continuous operation
Supply voltage:	100 V AC / 110 V AC 220 V AC 240 V AC
Mains frequency:	50 / 60 Hz
Power consumption:	max. 55 VA Scaler = 14 VA 27 - 30 kHz
Fuses:	Ø 5 x 20 mm, 250 V AC, time lag for 100 V AC / 110 V AC power supply: T 1 A for 220 V AC / 240 V AC power supply: T 500 mA L
Weight:	Approximately 3,6 kg
Dimensions:	Height: 120 mm Width: 220 mm Length: 240 mm
Operating conditions:	+10°C - +40°C, 30% - 75% relative humidity
Storage and transport conditions:	-10°C - +40°C, 10% - 95% relative humidity, 500 hPa - 1060 hPa air pressure

* For units sold in the USA or Canada.

Symbols

Front panel

"Water" (5) and "Power" (6) knob settings

 minimum

 maximum


Symbols on the rear panel

 Power switch "ON"

 Power switch "OFF"

 Caution! Read the operating instructions

 Applied part type B

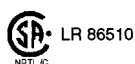CSA symbol marked "NRTL/C" for certified products
 (*applies to 100/110 V AC units only)
 CSA number (*applies to 100/110 V AC units only)

The fuse must only be replaced with the type indicated(*applies to 100/110 V AC units only)

CE symbol
 Refers to EU directive 93/42 EU
 Including EN 60601-1 and EN 60601-1-2

Marking on the side of the housing

Rocker switch

Marking on footswitch

IP21 Degree of protection against water

Text on the top of the housing

Rinse out with water after each use.
 (see booklet)

* For units sold in the USA or Canada.

Electromagnetic compatibility

Electromagnetic emissions

The Piezon Master 400 is intended for use in the electromagnetic environment specified below. The customer or the user of the Piezon Master 400 should assure that it is used in such an environment.

Emissions test	Compliance	Electromagnetic environment - guidance
RF emissions CISPR 11	Groupe 1	The Piezon Master 400 uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment. The Piezon Master 400 is suitable for use in all establishments, including domestic establishments and those directly connected to the public low-voltage power supply network that supplies buildings used for domestic purposes.
RF emissions CISPR 11	Class B	
Harmonic emissions IEC 61000-3-2	Class A	
Voltage fluctuations/Flicker emissions IEC 61000-3-3	Complies	

Electromagnetic immunity

The Piezon Master 400 is intended for use in the electromagnetic environment specified below. The customer or the user of the Piezon Master 400 should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Electrostatic discharge (ESD) IEC 61000-4-2	± 6 kV contact ± 8 kV air	± 6 kV contact ± 8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Electrical fast transient/burst IEC 61000-4-4	± 2 kV for power supply lines	± 2 kV for power supply lines	Mains power quality should be that of a typical commercial or hospital environment.
Surge IEC 61000-4-5	± 1 kV differential mode ± 2 kV common mode	± 1 kV differential mode ± 2 kV common mode	Mains power quality should be that of a typical commercial or hospital environment.
Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11	$<5\%$ U_T ($>95\%$ dip in U_T) for 0.5 cycle 40% U_T (60 % dip in U_T) for 5 cycles 70% U_T (30 % dip in U_T) for 25 cycles $<5\%$ U_T ($>95\%$ dip in U_T) for 5 sec	$<5\%$ U_T ($>95\%$ dip in U_T) for 0.5 cycle 40% U_T (60 % dip in U_T) for 5 cycles 70% U_T (30 % dip in U_T) for 25 cycles $<5\%$ U_T ($>95\%$ dip in U_T) for 5 sec	Mains power quality should be that of a typical commercial or hospital environment. If the user of the Piezon Master 400 requires continued operation during power mains interruptions, it is recommended that the Piezon Master 400 be powered from an uninterruptible power supply or a battery.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

NOTE: U_T is the a.c. mains voltage prior to application of the test level.

Electromagnetic immunity, RF portable equipment

The Piezon Master 400 is intended for use in the electromagnetic environment specified below. The customer or the user of the Piezon Master 400 should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment - guidance
Portable and mobile RF communications equipment should be used no closer to any part of the Piezon Master 400, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.			
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	3 Vrms	Recommended separation distance $d = 1.2 \sqrt{P}$

NOTE 1: At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which the Piezon Master 400 is used exceeds the applicable RF compliance level above, the Piezon Master 400 should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as re-orienting or relocating the Piezon Master 400.

b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

Recommended separation distances

The Piezon Master 400 is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of the Piezon Master 400 can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and the Piezon Master 400 as recommended below, according to the maximum output power of the communications equipment.

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter m		
	150 kHz to 80 MHz $d = 1.2 \sqrt{P}$	80 MHz to 800 MHz $d = 0.2 \sqrt{P}$	800 MHz to 2,5 GHz $d = 0.4 \sqrt{P}$
0.01	0.12 m	0.12 m	0.23 m
0.1	0.38 m	0.38 m	0.73 m
1	1.2 m	1.2 m	2.3 m
10	3.8 m	3.8 m	7.3 m
100	12 m	12 m	23 m
For transmitters rated at a maximum output power not listed above, the recommended separation distance d in metres (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.			

NOTE 1: At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2: These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Cables length

Cables and accessories	Maximum length	Complies with
Handpiece cord	1.7 m	RF emissions, CISPR 11, Class B/Group 1
Footswitch cable	3 M	Harmonic emissions, IEC 61000-3-2 Voltage fluctuations/flicker emission, IEC 61000-3-3 Electrostatic discharge (ESD), IEC 61000-4-2 Electric fast transient/burst, IEC 61000-4-4 Surge, IEC 61000-4-5 Voltage dips, short interruptions and voltage variations on power supply input lines IEC 61000-4-11 Power frequency (50/60 Hz) magnetic field IEC 61000-4-8 Conducted RF IEC 61000-4-6 Radiated RF, IEC 61000-4-3

Troubleshooting

Fault	Possible cause	Remedy
The green equipment-on lamp (3) does not illuminate when the unit is switched on	No power supply, power switch (2) OFF	Switch the power switch (2) ON
	Mains plug not plugged in or defective	Check the mains plug
	Mains socket not live	Check the power supply in the practice
	The green lamp is defective	Send the unit to the EMS service department
	Defective fuse	Send the unit to the EMS service department
	Electrical fault	Send the unit to the EMS service department
No ultrasonic oscillations when the footswitch (18) is activated	Power switch (2) switched OFF	Switch power switch (2) ON
	Use the first footswitch setting with the rocker switch set to IRRIGATION/SCALING	Press on the footswitch firmly
	The handpiece hose (8) is in the handpiece holder (9)	Remove the handpiece (8) from the handpiece holder (9)
	The instrument or file holder have not been fixed in place properly	Check that the instrument or file holder has been fixed in place properly
	The handpiece of the system has not been fixed in place properly	Check that the handpiece has been fixed in place properly
	There is liquid in the handpiece hose (8) coupling	Dry the O-ring with compressed air and check the handpiece plug
	The instrument or file holder is damaged	Order a new instrument or file holder
	The handpiece hose (8) or handpiece is defective	Send the unit to the EMS service department
	Defective electronics	Send the unit to the EMS service department

Fault	Possible cause	Remedy
No irrigating liquid flows when the footswitch (18) is activated	Power switch (2) switched OFF	Switch power switch (2) ON
	"Water" knob (5) is not set high enough	Set "Water" knob (5) higher
	The bottle (10) is not in place or is empty	Fill the bottle (10) and fit it
	The handpiece hose (8) is in the handpiece holder (9)	Remove the handpiece hose (8) from the handpiece holder (9)
	Use the second footswitch setting with the rocker switch set to ENDO	Only press the footswitch gently
	A dry work handpiece has been fitted	Change the handpiece
	The instrument is blocked	Clean the instrument with compressed air
	The sections where the liquid passes through the unit, handpiece hose or handpiece are blocked	Send the unit and handpiece to the EMS service department
	The handpiece hose (8) or handpiece are defective	Send the unit to the EMS service department
Liquid leaks from the irrigating liquid bottle	The red washer (13) has been fitted to the cap incorrectly, has been damaged or is missing	Check the red washer (13)
	The screw-cap of the bottle (11) has been damaged	Replace the screw-cap
	The O-ring on the small nipple in the bottle holder (15) has been damaged	Replace the O-rings
	The bayonet face plate has been damaged	Replace the bayonet face plate
	The irrigating liquid bottle has been damaged	Replace the bottle
Liquid leaks from the handpiece hose coupling	The O-ring on the handpiece plug is defective	Replace the O-ring
Liquid leaks from the unit	A component where the liquid passes through in the unit has been damaged	Send the unit to the EMS service department

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EMS ELECTRO MEDICAL SYSTEMS VERTRIEBS GmbH, Schatzbogen 86, D-81829 München